




## 4-Methoxy-phencyclidine

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.

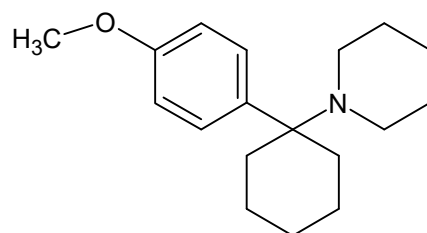

### 1. GENERAL INFORMATION

**IUPAC Name:** 1-[1-(4-methoxyphenyl)cyclohexyl]piperidine

**CAS#:**

**Synonyms:** 4-MeO PCP, 4-MeO Phencyclidine, 4-Methoxy PCP, Methoxydine

**Source:** DEA Reference Material Collection

**Appearance:** White powder

**UV<sub>max</sub>(nm):** Not determined

### 2. CHEMICAL AND PHYSICAL DATA

#### 2.1 CHEMICAL DATA

Form	Chemical Formula	Molecular Weight	Melting Point (°C)
Base	C <sub>18</sub> H <sub>27</sub> NO	273	Not Determined
HCl	C <sub>18</sub> H <sub>27</sub> NO HCl	310	184.6



## 4-Methoxy-phencyclidine

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.




### 3. QUALITATIVE DATA

#### 3.1 NUCLEAR MAGNETIC RESONANCE

**Sample Preparation:** Dilute analyte to ~7 mg/mL in DMSO- $d_6$  containing TMS for 0 ppm reference and maleic acid as quantitative internal standard.

**Instrument:** 400 MHz NMR spectrometer  
**Parameters:** Spectral width: at least containing -3 ppm through 13 ppm  
Pulse angle: 90°  
Delay between pulses: 45 seconds

$^1\text{H}$ NMR: 4-Methoxy-phencyclidine HCl; Lot# 0444879-31; DMSO- $d_6$ ; 400MHz

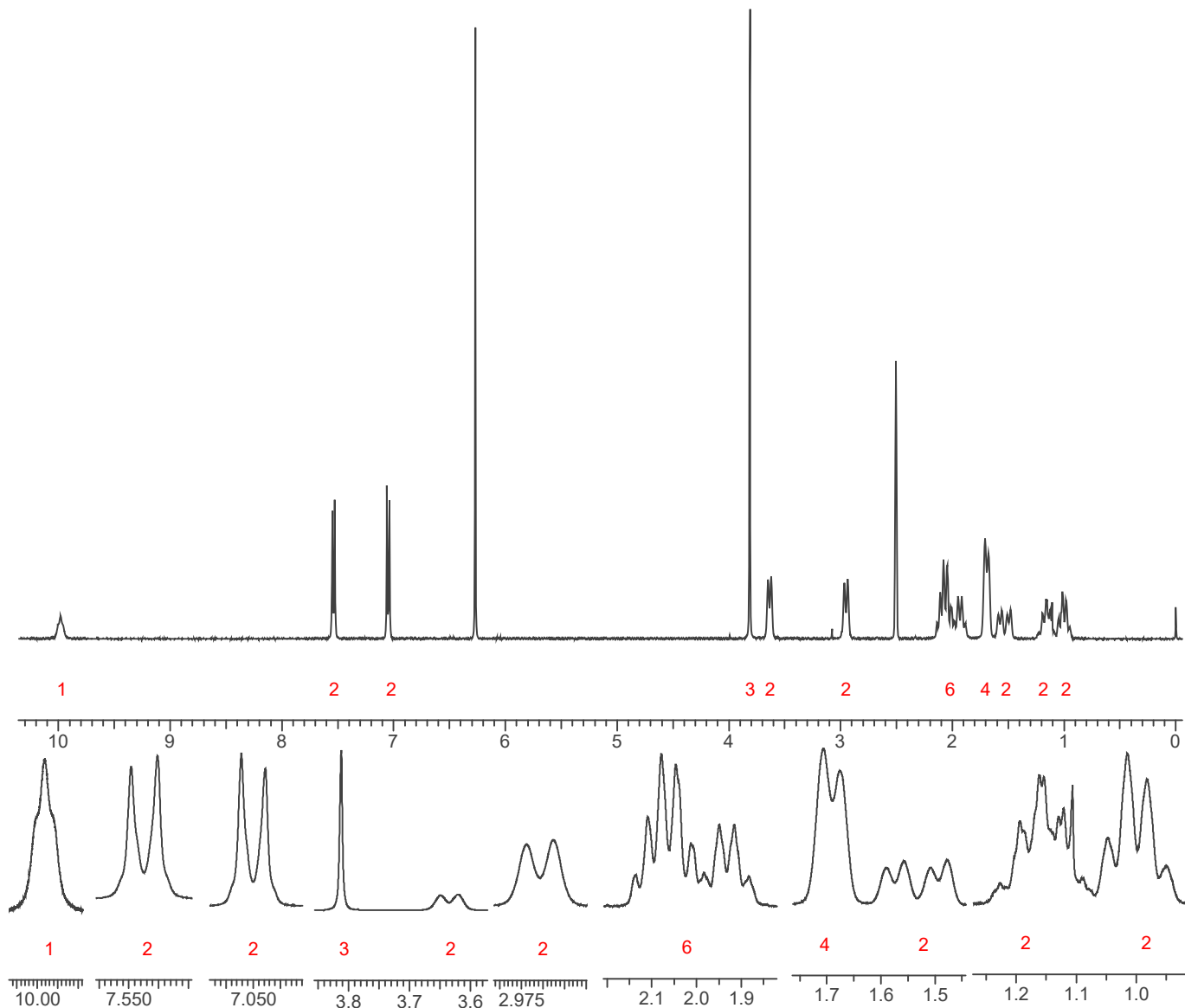



## 4-Methoxy-phencyclidine

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.




### 3.2 GAS CHROMATOGRAPHY/MASS SPECTROMETRY

*Sample Preparation:* Dilute analyte ~4 mg/mL in MeOH

**Instrument:** Agilent gas chromatograph operated in split mode with MS detector  
**Column:** HP-5 MS (or equivalent); 30m x 0.25 mm x 0.25  $\mu$ m  
**Carrier Gas:** Helium at 1.5 mL/min  
**Temperatures:** Injector: 280°C                      MSD transfer line: 280°C  
MS Source: 250°C                      MS Quad: 150°C  
Oven program:  
    1) 100°C initial temperature for 1.0 min  
    2) Ramp to 300°C at 12 °C/min  
    3) Hold final temperature for 9.0 min  
**Injection Parameters:** Split Ratio = 25:1, 1  $\mu$ L injected  
**MS Parameters:** Mass scan range: 30-550 amu                      Threshold: 150  
Tune file: stune.u                      Acquisition mode: scan  
**Retention Time:** 13.840 min

EI Mass Spectrum: 4-Methoxy-phencyclidine HCl; Lot# 0444879-31

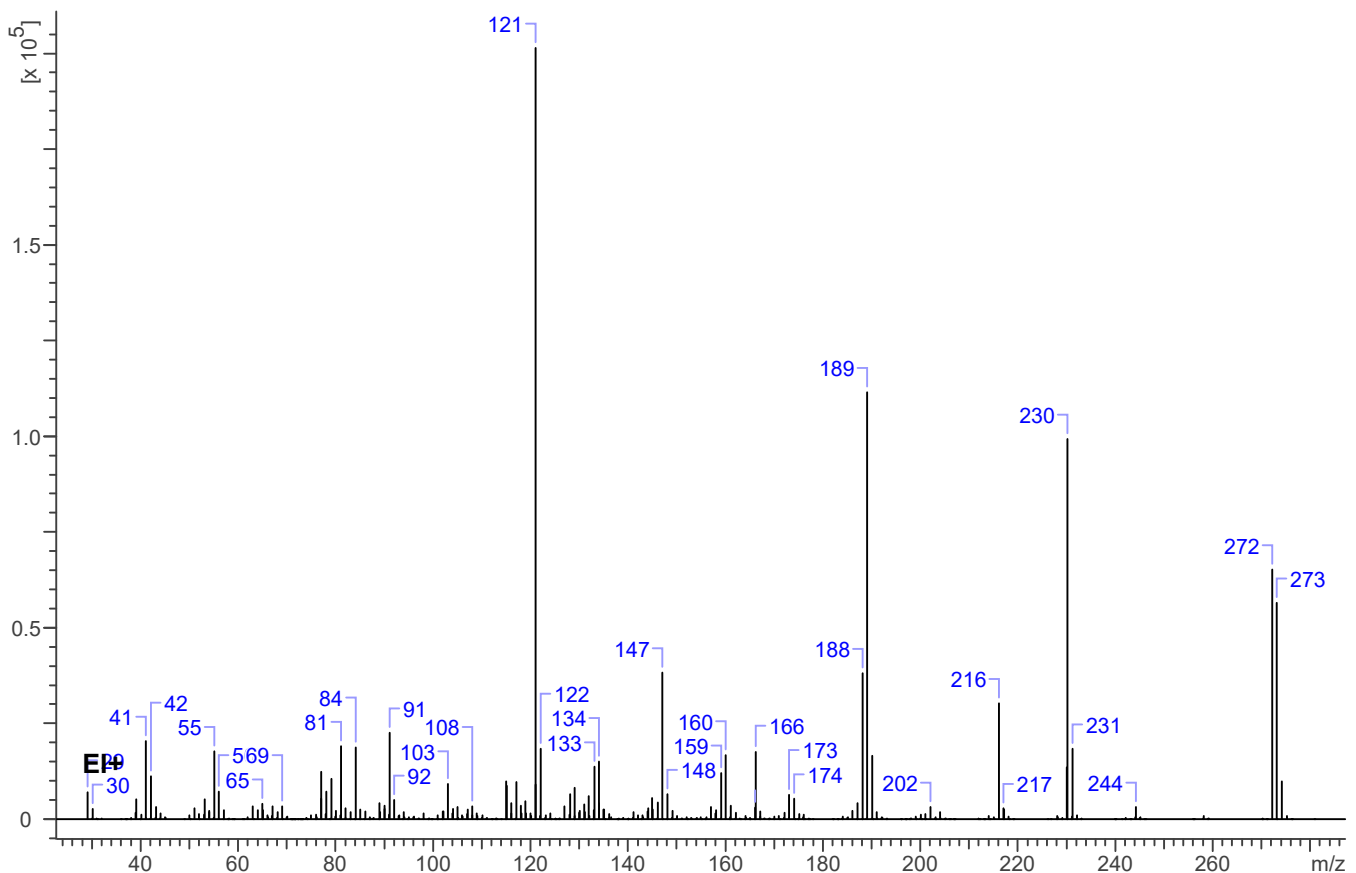



# 4-Methoxy-phencyclidine

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.




## 3.3 INFRARED SPECTROSCOPY (FTIR)

**Instrument:** FTIR with diamond ATR attachment (1 bounce)

**Scan Parameters:**  
Number of scans: 32  
Number of background scans: 32  
Resolution: 4 cm<sup>-1</sup>  
Sample gain: 8  
Aperture: 150

FTIR ATR (Diamond 1 Bounce): 4-Methoxy-phencyclidine HCl; Lot# 0444879-31

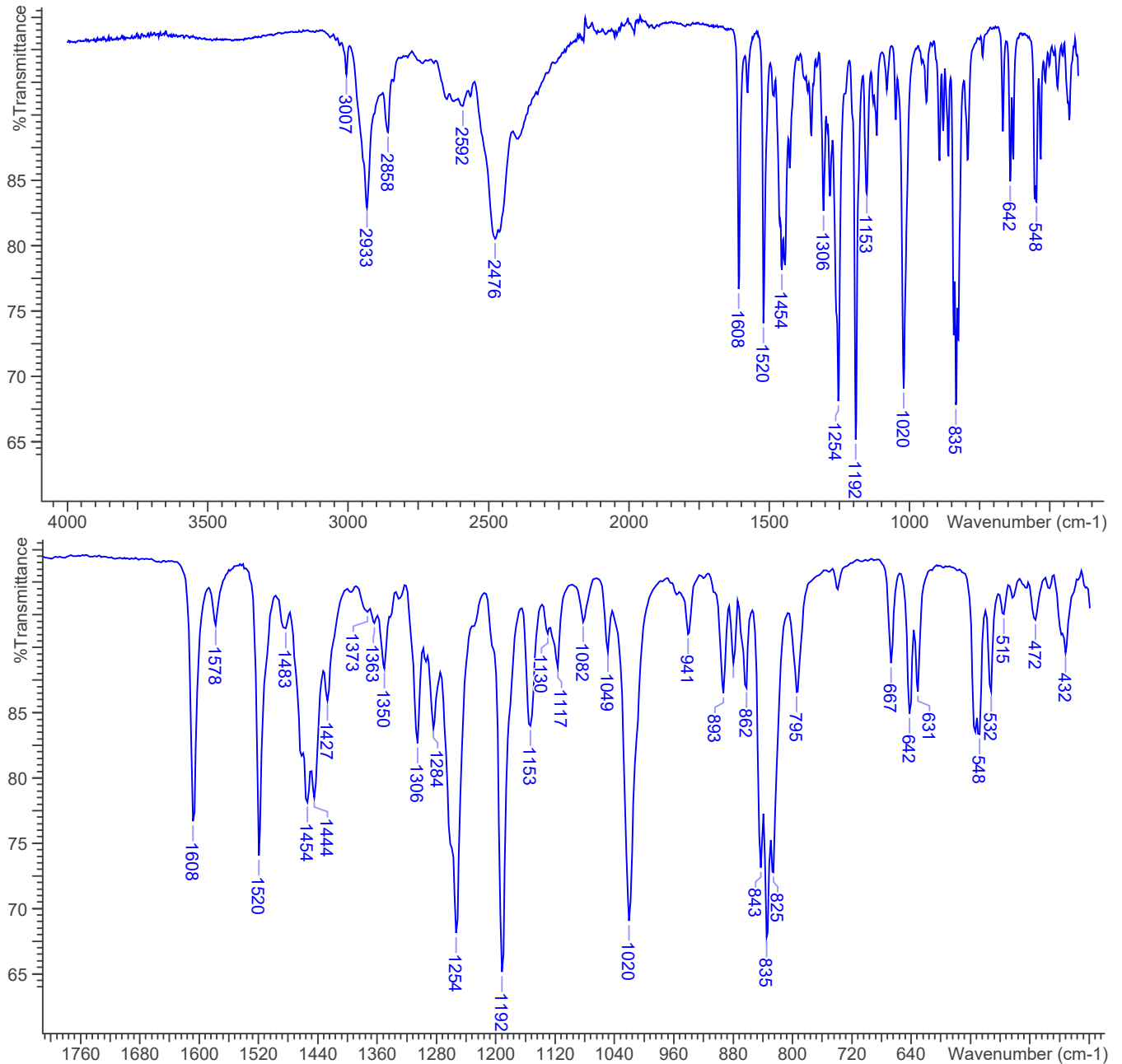



## 4-Methoxy-phencyclidine

*The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.*



### 4. ADDITIONAL RESOURCES

[\*Wikipedia\*](#)

[\*Forendex\*](#)