


AH-7563 hydrochloride

The Krstenansky lab at the KGI School of Pharmacy and Health Sciences generated this monograph using synthesized material

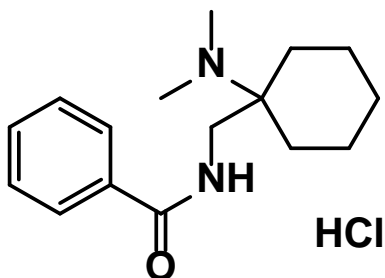

1. GENERAL INFORMATION

IUPAC Name: 1-(benzamidomethyl)-cyclohexyldimethylamine; hydrochloride

CAS#: 41804-98-2; 63886-94-2 (base)

Synonyms: A13

Source: Synthesized Material Lot# JLK010-027-A13

Appearance: White Crystals (HCl)

UV_{max} (nm): Not Determined

2. CHEMICAL AND PHYSICAL DATA

2.1 CHEMICAL DATA

Form	Chemical Formula	Molecular Weight	Melting Point (°C)
HCl	C ₁₆ H ₂₄ N ₂ O·HCl	296.84	232.6 ± 0.06
base	C ₁₆ H ₂₄ N ₂ O	260.37	Not determined

3. QUALITATIVE DATA

3.1 NUCLEAR MAGNETIC RESONANCE

Sample Preparation: Dilute analyte to ~5 mg/mL in deuterated chloroform:methanol (CDCl₃:CD₃OD, 1:5) + TMS.



Instrument: 400 MHz NMR spectrometer

Parameters: Spectral width: 6410.3 Hz containing -3 ppm through 13 ppm

Pulse angle: 90°

Delay between pulses: 30 seconds

¹H NMR: A13 HCl; Lot JLK010-027-A13; CDCl₃:CD₃OD (1:5) + TMS; 400 MHz

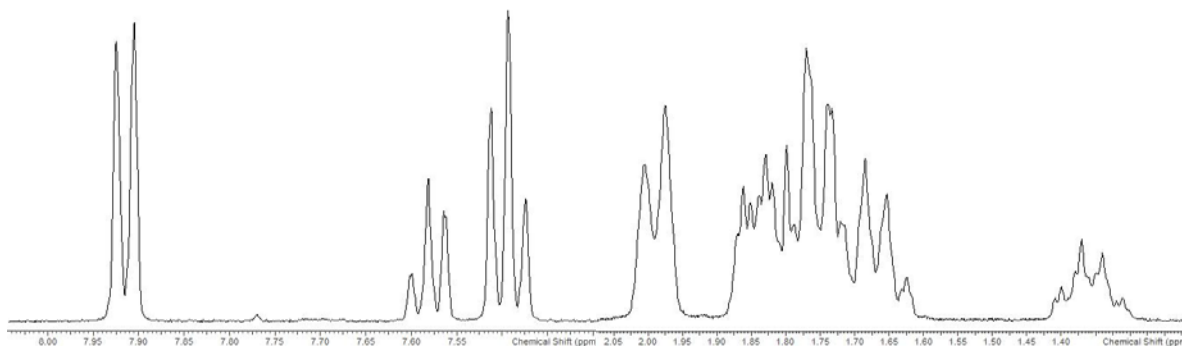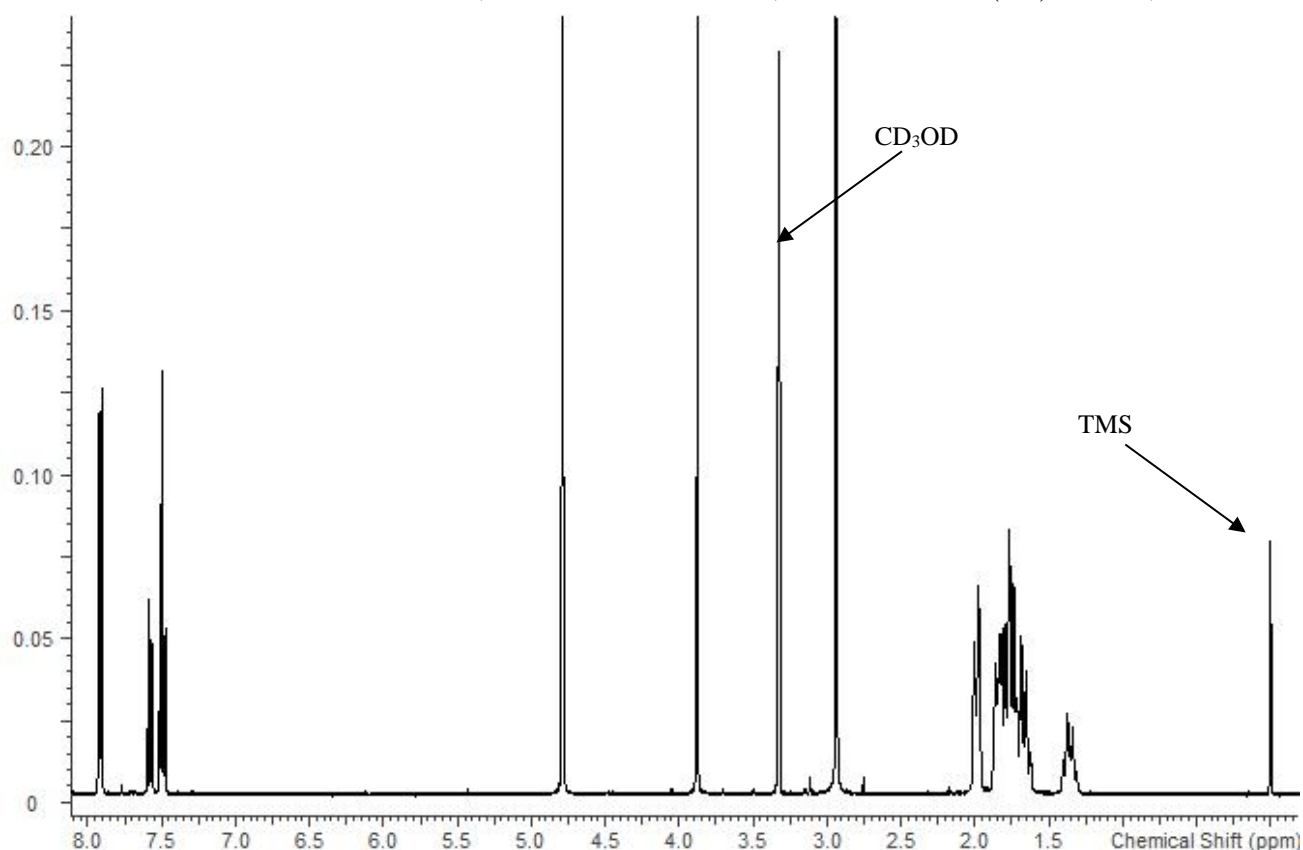


3.2 GAS CHROMATOGRAPHY/MASS SPECTROMETRY


Sample Preparation: Dilute analyte ~ 1 mg/mL in methanol

Instrument:	Shimadzu gas chromatograph operated in split mode with MS detector
Column:	Rtx5MS (a DB-5 equivalent); 30m x 0.25 mm x 0.25 μm
Carrier Gas:	Helium at 1 mL/min
Temperatures:	Injector: 280°C MSD transfer line: 280°C MS Source: 200°C Oven program: 1) 90°C initial temperature for 2.0 min 2) Ramp to 300°C at 14°C/min 3) Hold final temperature for 10.0 min
Injection Parameters:	Split Ratio = 1:15, 1 μL injected
MS Parameters:	Mass scan range: 34-550 amu Threshold: 100 Tune file: 050218_Tune.qgt Acquisition mode: scan
Retention Time:	14.82 min

EI Mass Spectrum: A13 HCl; Lot JLK010-027-A13

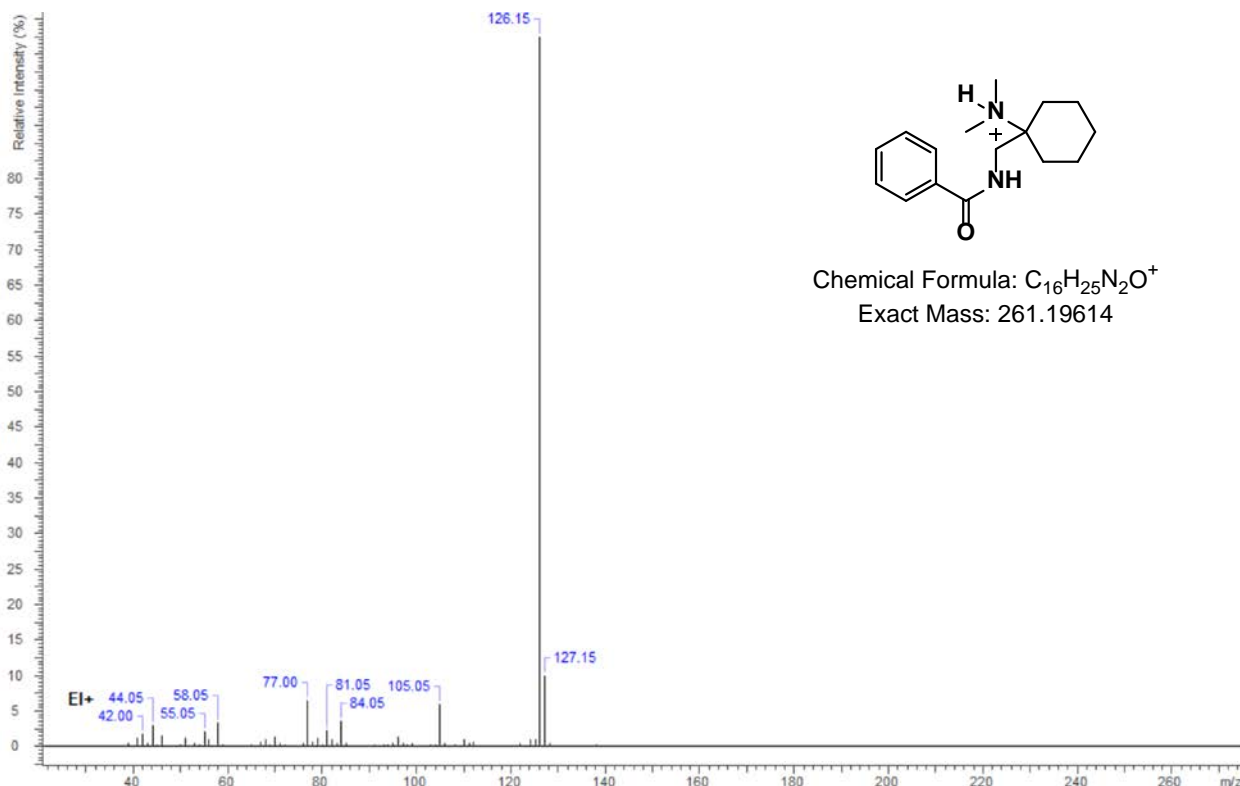

Zoomed view (126.15 is 100% relative intensity and is truncated in this view)

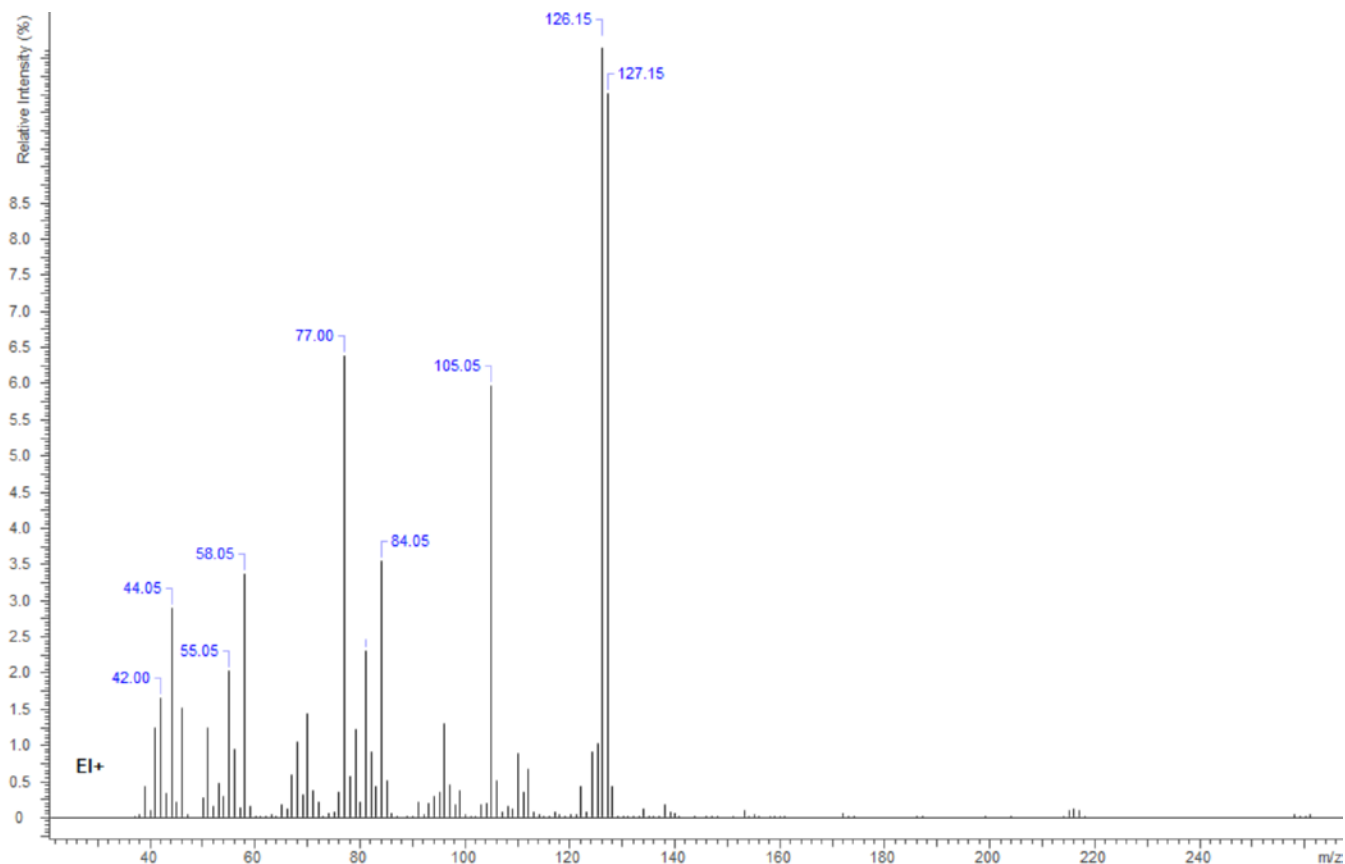


3.3 INFRARED SPECTROSCOPY (FTIR)

Instrument: FTIR with ZnSe ATR attachment (1 bounce)


Scan Parameters: Number of scans: 4
 Number of background scans: 4
 Resolution: 4 cm⁻¹
 Sample gain: 8
 Aperture: 150

FTIR ATR (ZnSe, 1 Bounce): A13 HCl; Lot JLK010-027-A13

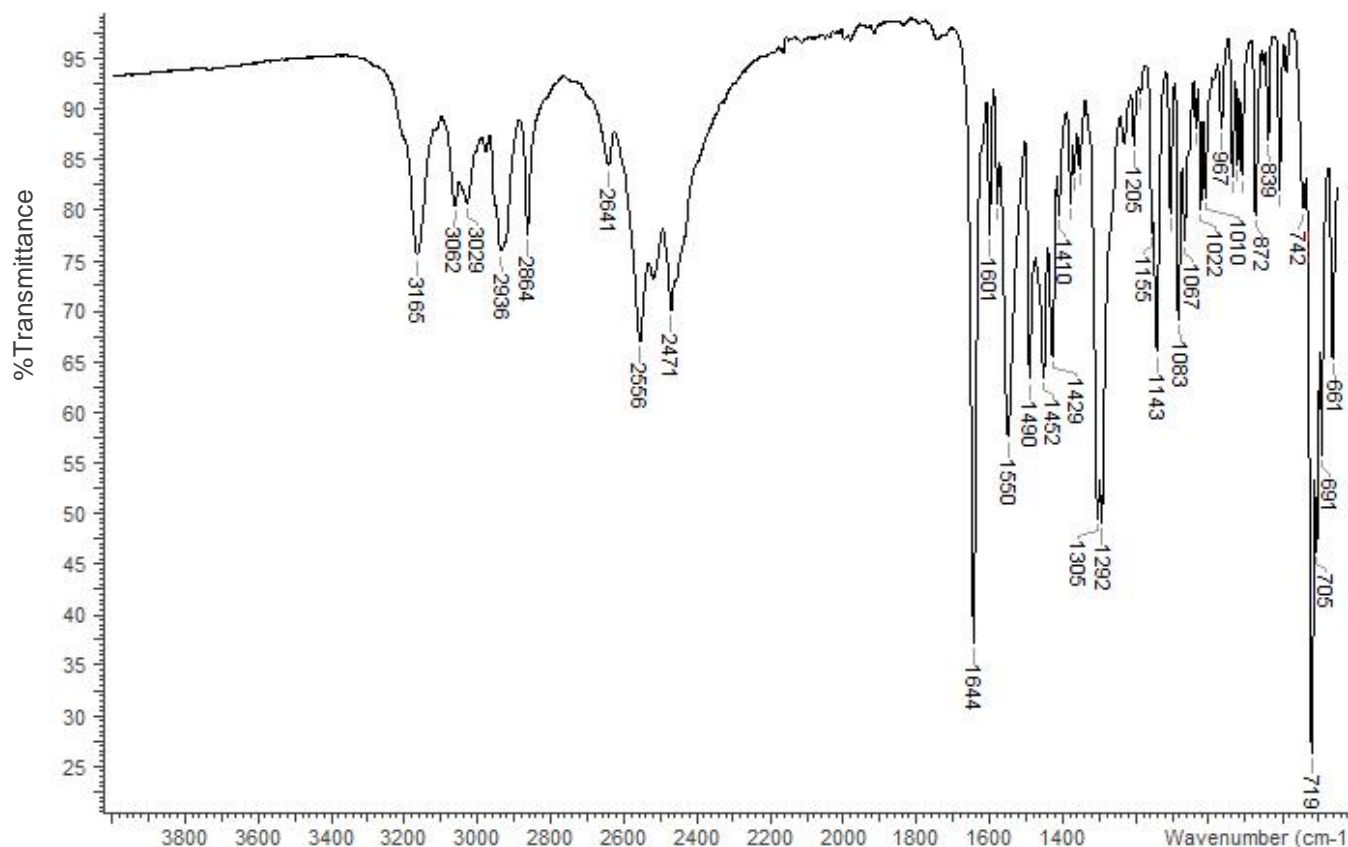

AH-7563 hydrochloride


The Krstenansky lab at the KGI School of Pharmacy and Health Sciences generated this monograph using synthesized material

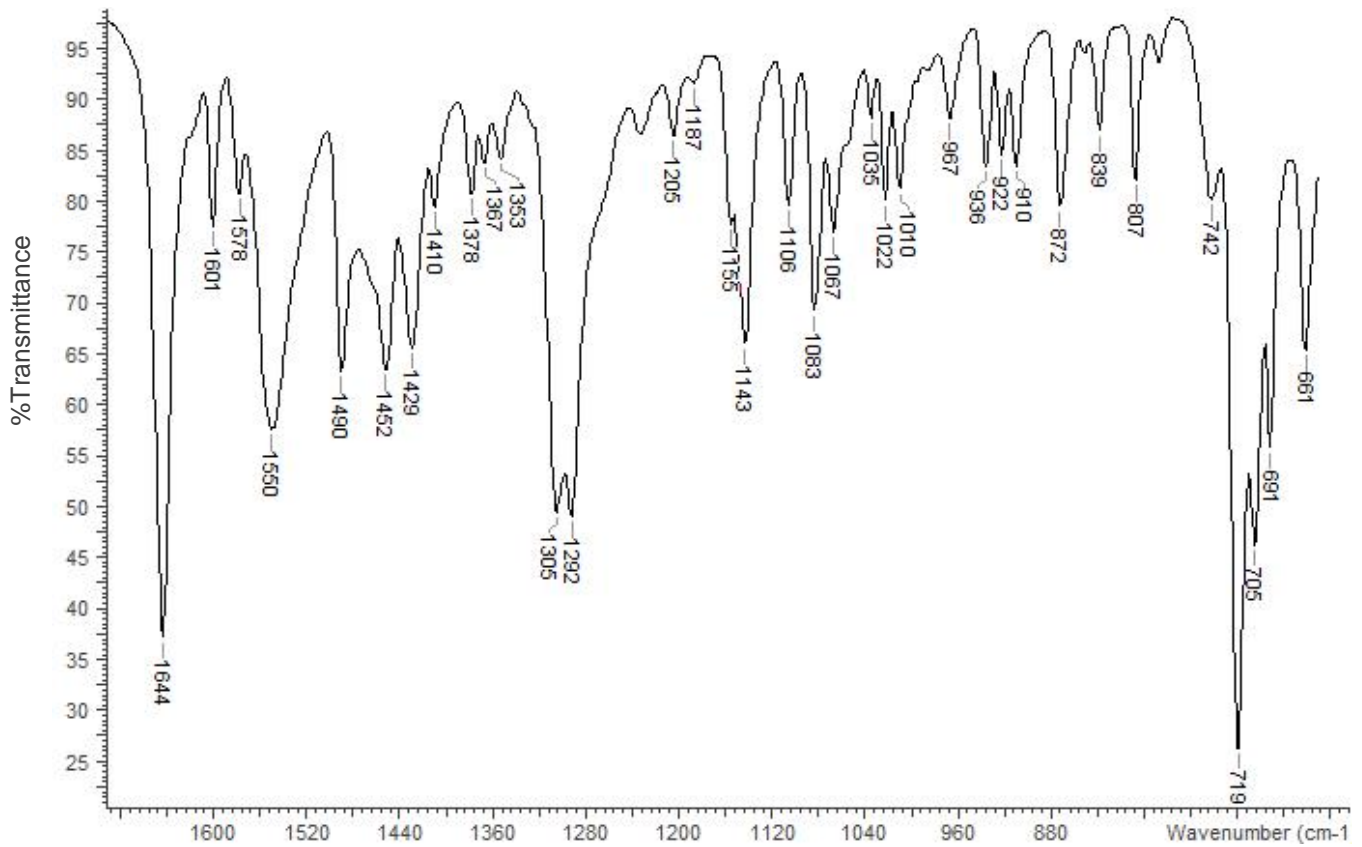

3.4 RAMAN SPECTROSCOPY

Instrument: Rigaku Progeny 1064
Scan Parameters: Power (mW): 350
Exposure (ms): 1000
Averages: 30
Threshold: 0.80

Raman (1064 nm): A13 HCl; Lot JLK010-027-A13

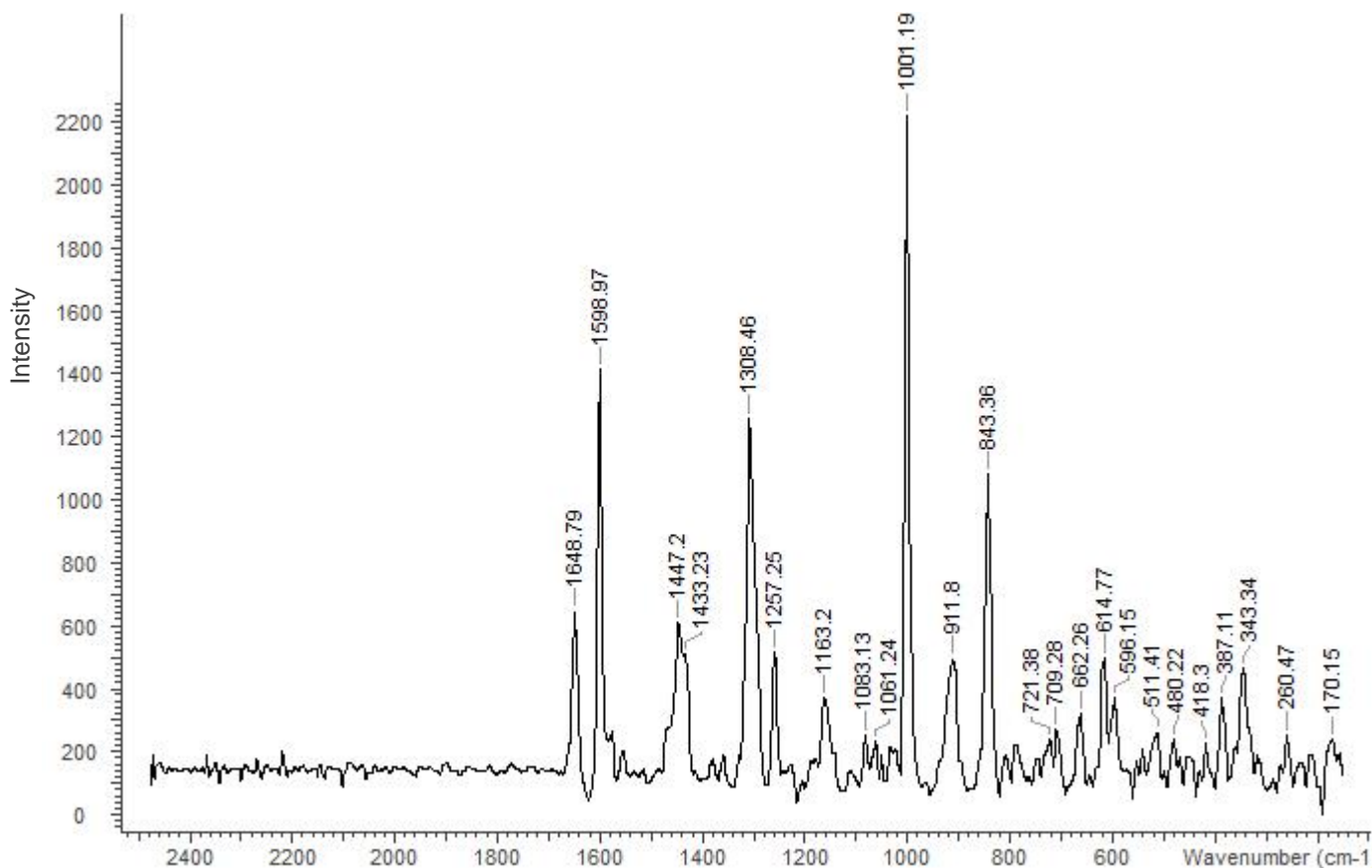

4. ADDITIONAL RESOURCES

1-(3,4-DICHLOROBENZAMIDOMETHYL)CYCLOHEXYLDIMETHYLAMINE

Norman James Harper and George Bryan Austin Veitch

US Patent 3,975,443 Aug. 17, 1976 Example 4

1-(3,4-Dichlorobenzamidomethyl)cyclohexyldimethylamine and related compounds as potential analgesics

N. J. Harper, G. B. A. Veitch, and D. G. Wibberley

Journal of Medicinal Chemistry 1974 17 (11), 1188-1193

DOI: 10.1021/jm00257a012

Compound AH-7563

Tom Hsu, Jayapal Reddy Mallareddy, Kayla Yoshida, Vincent Bustamante, Tim Lee, John L. Krstenansky, Alexander C. Zambon, Synthesis and pharmacological characterization of ethylenediamine synthetic opioids in human μ -opiate receptor 1 (OPRM1) expressing cells. *Pharmacol. Research & Perspectives* 7: e00511 (2019) doi: 10.1002/prp2.511

5. ACKNOWLEDGEMENT

These data are from a project supported by Award No. 2016-R2-CX-0059, awarded by the National Institute of Justice, Office of Justice Programs, U.S. Department of Justice. The opinions, findings, and conclusions or recommendations expressed in this publication are those of the authors and do not necessarily reflect those of the Department of Justice. We also thank Rigaku Corporation for the loan of the Progeny 1064 Raman instrument.