




# Flurazepam

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.

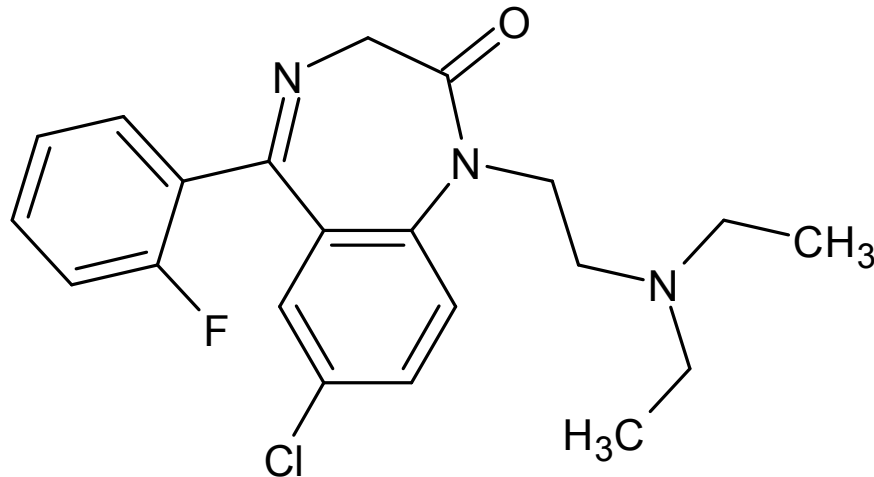

## 1. GENERAL INFORMATION

<b>IUPAC Name:</b>	7-chloro-1-[2-(diethylamino)ethyl]-5-(2-fluorophenyl)-1,3-dihydro-2H-1,4-benzodiazepin-2-one
<b>CAS#:</b>	17617-23-1
<b>Synonyms:</b>	Dalmane, Dalmadorm, Benozil
<b>Source:</b>	DEA Reference Material Collection
<b>Appearance:</b>	White powder (diHCl)
<b>UV<sub>max</sub>(nm):</b>	Not Determined

## 2. CHEMICAL AND PHYSICAL DATA

### 2.1 CHEMICAL DATA

Form	Chemical Formula	Molecular Weight	Melting Point (°C)
Base	C <sub>21</sub> H <sub>23</sub> ClFN <sub>3</sub> O	387	Not Determined
diHCl	C <sub>21</sub> H <sub>23</sub> ClFN <sub>3</sub> O · 2HCl	460	213.2



# Flurazepam

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.




## 3. QUALITATIVE DATA

### 3.1 NUCLEAR MAGNETIC RESONANCE

**Sample Preparation:** Dilute analyte to ~10 mg/mL in CD<sub>3</sub>OD containing TMS for 0 ppm reference and dimethylsulfone as quantitative internal standard.

**Instrument:** 400 MHz NMR spectrometer  
**Parameters:** Spectral width: at least containing -3 ppm through 13 ppm  
Pulse angle: 90°  
Delay between pulses: 45 seconds

<sup>1</sup>H NMR: Flurazepam diHCl Lot # SLBG2349V; CD<sub>3</sub>OD; 400MHz

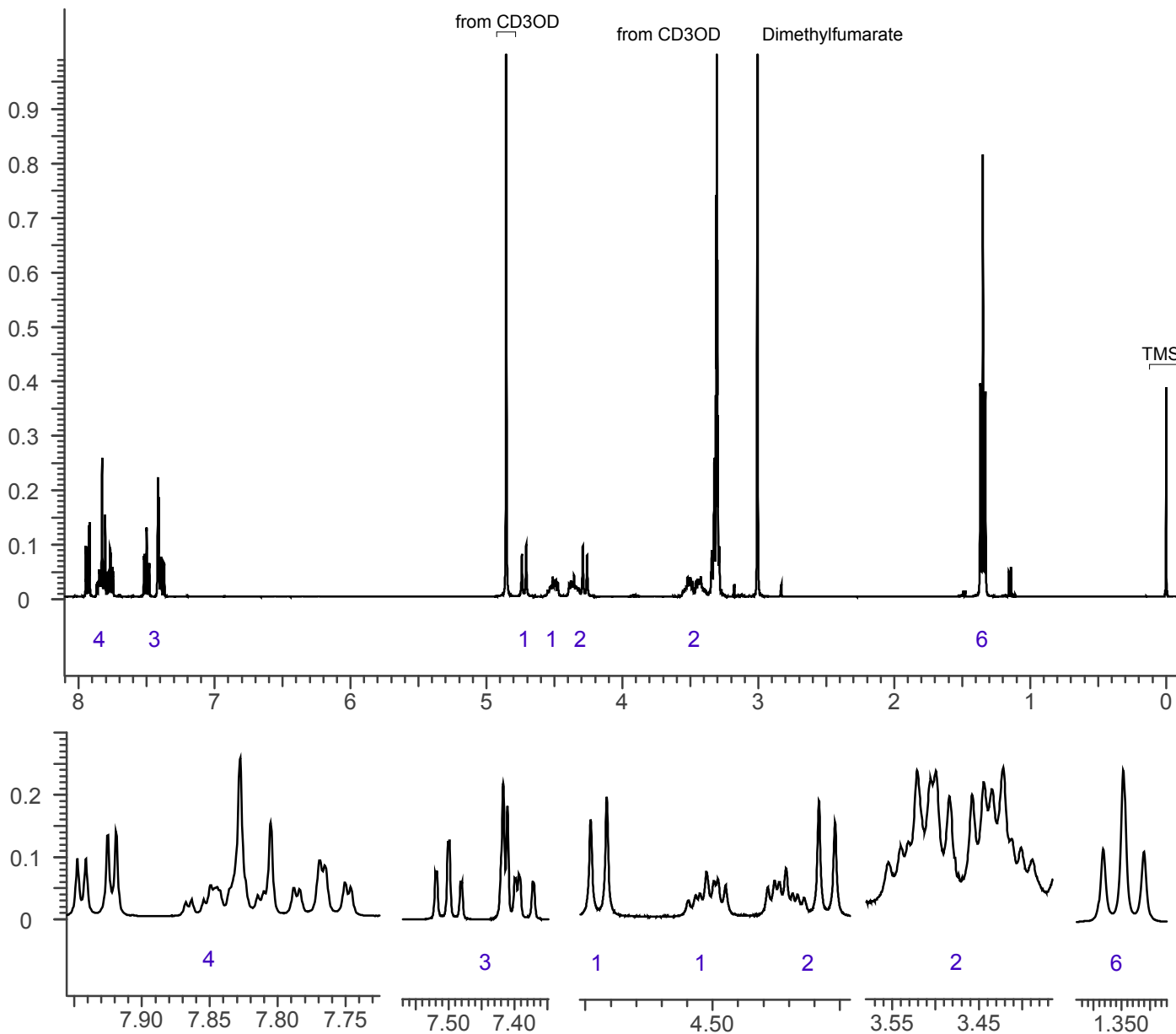



# Flurazepam

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.



## 3.2 GAS CHROMATOGRAPHY/MASS SPECTROMETRY


*Sample Preparation:* Dilute analyte ~4 mg/mL base extracted into chloroform.

**Instrument:** Agilent gas chromatograph operated in split mode with MS detector  
**Column:** DB-1 MS (or equivalent); 30m x 0.25 mm x 0.25  $\mu$ m  
**Carrier Gas:** Helium at 1 mL/min  
**Temperatures:** Injector: 280°C  
MSD transfer line: 280°C  
MS Source: 230°C  
MS Quad: 150°C  
Oven program:  
1) 100°C initial temperature for 1.0 min  
2) Ramp to 300°C at 12 °C/min  
3) Hold final temperature for 9.0 min

**Injection Parameters:** Split Ratio = 20:1, 1  $\mu$ L injected  
**MS Parameters:** Mass scan range: 30-550 amu  
Threshold: 100  
Tune file: stune.u  
Acquisition mode: scan

**Retention Time:** 17.482 min

EI Mass Spectrum: Flurazepam diHCl Lot # SLBG2349V

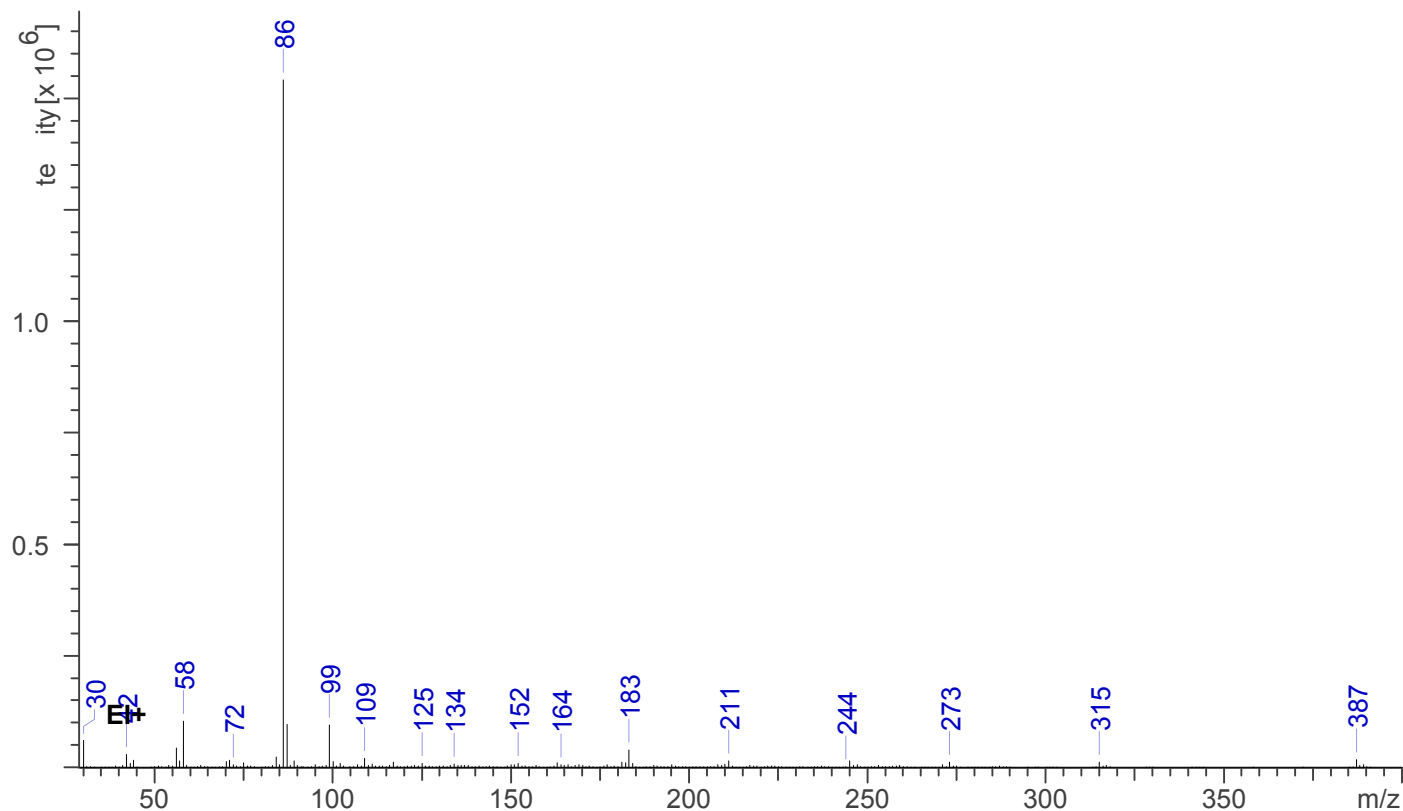



# Flurazepam




The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.

## 3.3 INFRARED SPECTROSCOPY (FTIR)

**Instrument:** FTIR with diamond ATR attachment (3 bounce)

**Scan Parameters:**  
Number of scans: 32  
Number of background scans: 32  
Resolution: 4 cm<sup>-1</sup>  
Sample gain: 8  
Aperture: 150

FTIR ATR (Diamond, 3 Bounce): Flurazepam diHCl Lot # SLBG2349V

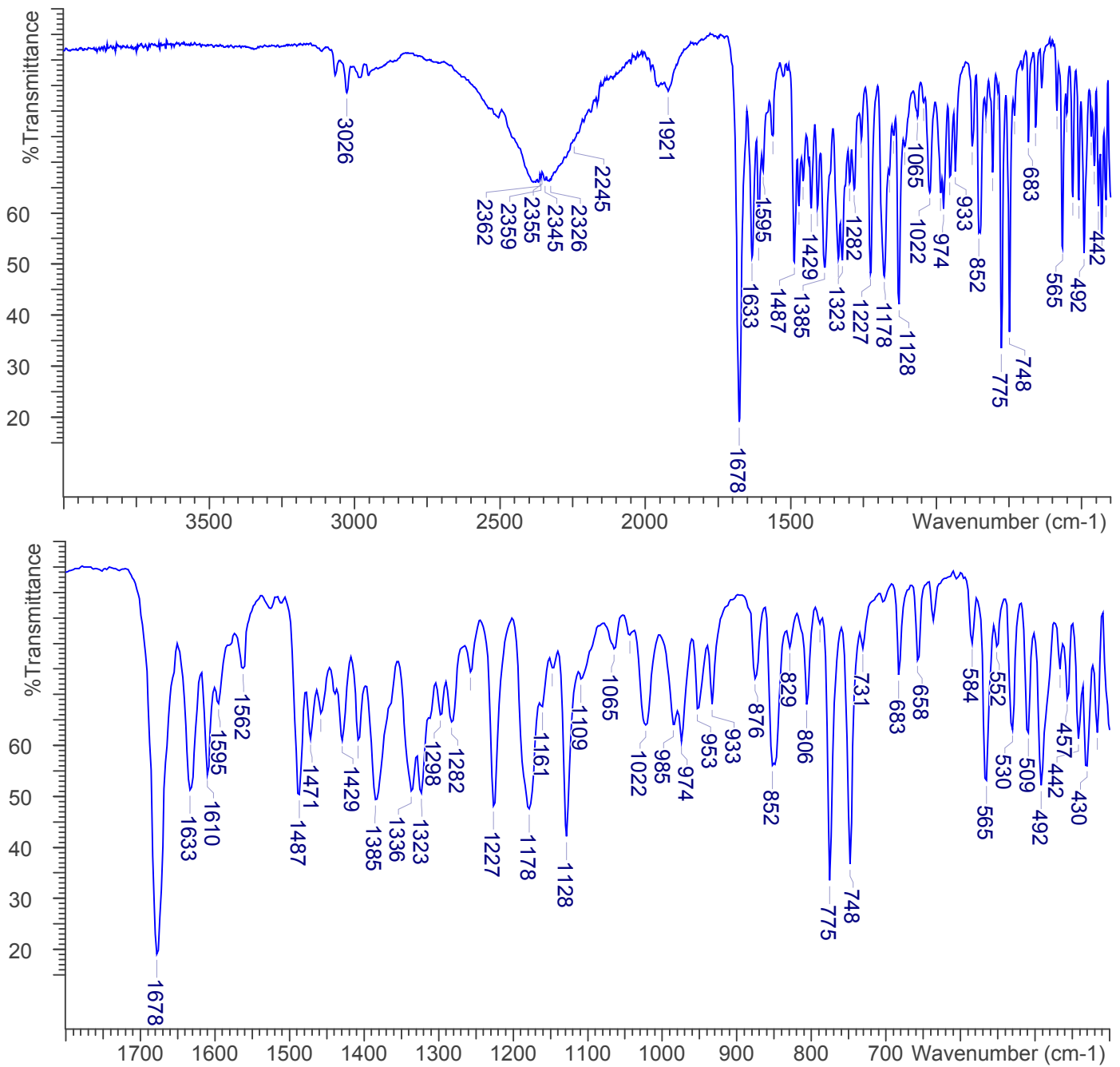



## Flurazepam

*The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.*



### **4. ADDITIONAL RESOURCES**

[Forendex](#)

[Wikipedia](#)