

1. GENERAL INFORMATION

| | |
|-------------------------------|-----------------------------------|
| IUPAC Name: | N-ethyl-1-phenylpropan-2-amine |
| CFR: | Schedule I |
| CAS#: | 457-87-4 (base), 16105-78-5 (HCl) |
| Synonyms: | Ethylamphetamine, Etilamfetamine |
| Source: | DEA Reference Material Collection |
| Appearance: | White powder (HCl) |
| Kovat's Index: | Pending |
| UV_{max} (nm): | Not Determined |

2. CHEMICAL AND PHYSICAL DATA

2.1 CHEMICAL DATA

| Form | Chemical Formula | Molecular Weight | Melting Point (°C) |
|------|---|------------------|--------------------|
| Base | C ₁₁ H ₁₇ N | 163 | Not Determined |
| HCl | C ₁₁ H ₁₇ N · HCl | 199 | 146.7 |

3. ADDITIONAL RESOURCES

[Forendex](#)

[Wikipedia](#)

4. QUALITATIVE DATA

4.1 NUCLEAR MAGNETIC RESONANCE

Method NMR D₂O

Sample Preparation: Dilute analyte to ~5 mg/mL in D₂O containing TSP for 0 ppm reference and maleic acid as quantitative internal standard.


Instrument: 400 MHz NMR spectrometer

Parameters: Spectral width: at least containing -3 ppm through 13 ppm

Pulse angle: 90°

Delay between pulses: 45 seconds

¹H NMR: N-ethylamphetamine HCl; Lot 2TDM-164-01; D₂O; 400 MHz

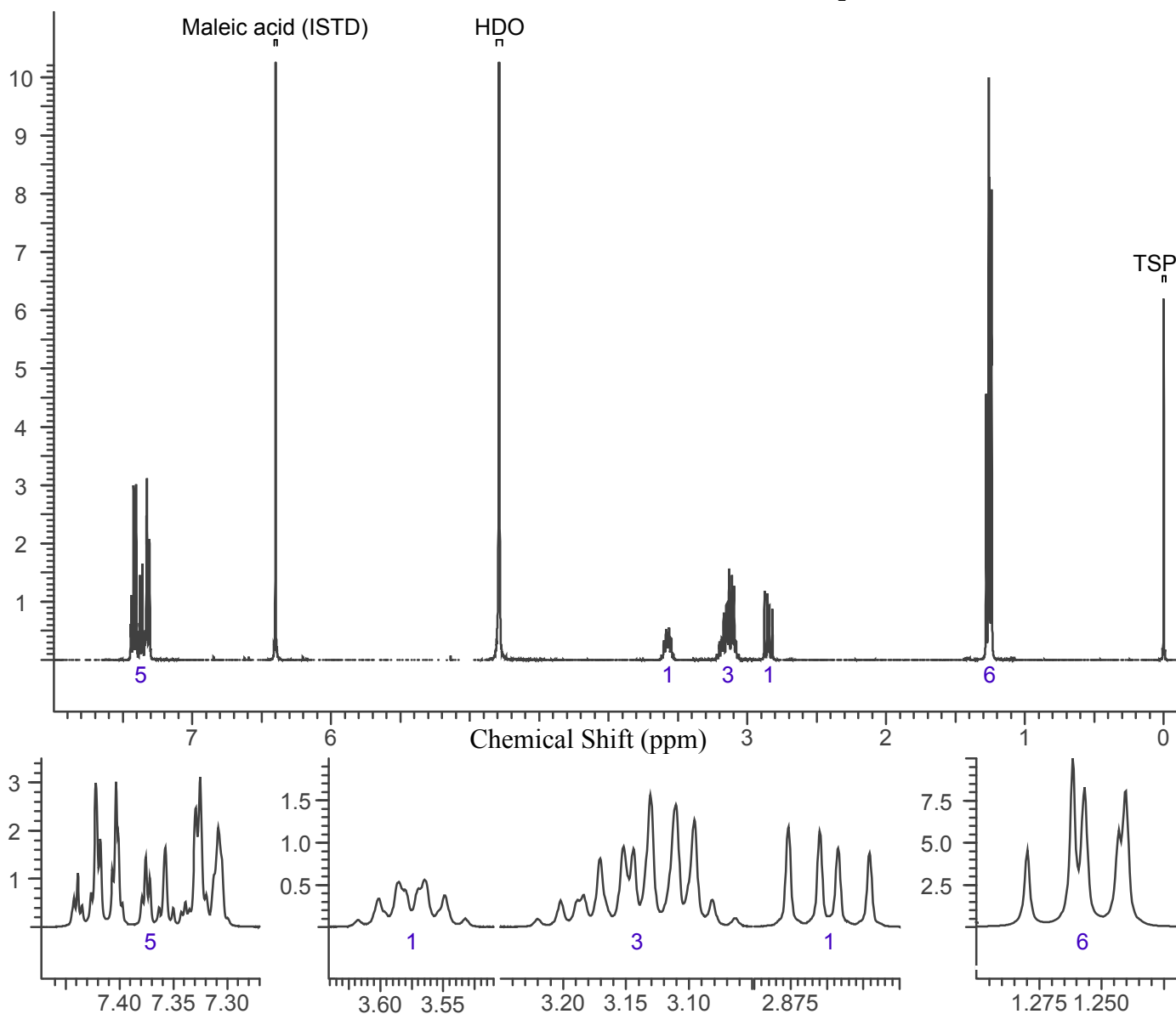

4.2 Gas Chromatography/Mass Spectrometry

Sample Preparation: Dilute analyte ~ 1 mg/mL base extracted into chloroform.

Instrument: Agilent gas chromatograph operated in split mode with MS detector

Column: DB-1 MS (or equivalent); 30m x 0.25 mm x 0.25 μ m

Carrier Gas: Helium at 1 mL/min

Temperatures: Injector: 280°C

MSD transfer line: 280°C

MS Source: 230°C

MS Quad: 150°C

Oven program:

1) 100°C initial temperature for 1.0 min

2) Ramp to 300°C at 12 °C/min

3) Hold final temperature for 9.0 min

Injection Parameters: Split Ratio = 20:1, 1 μ L injected

MS Parameters: Mass scan range: 30-550 amu


Threshold: 90

Tune file: stune.u

Acquisition mode: scan

Retention Time: 5.224 min

EI Mass Spectrum: N-ethylamphetamine HCl; Lot 2TDM-164-01

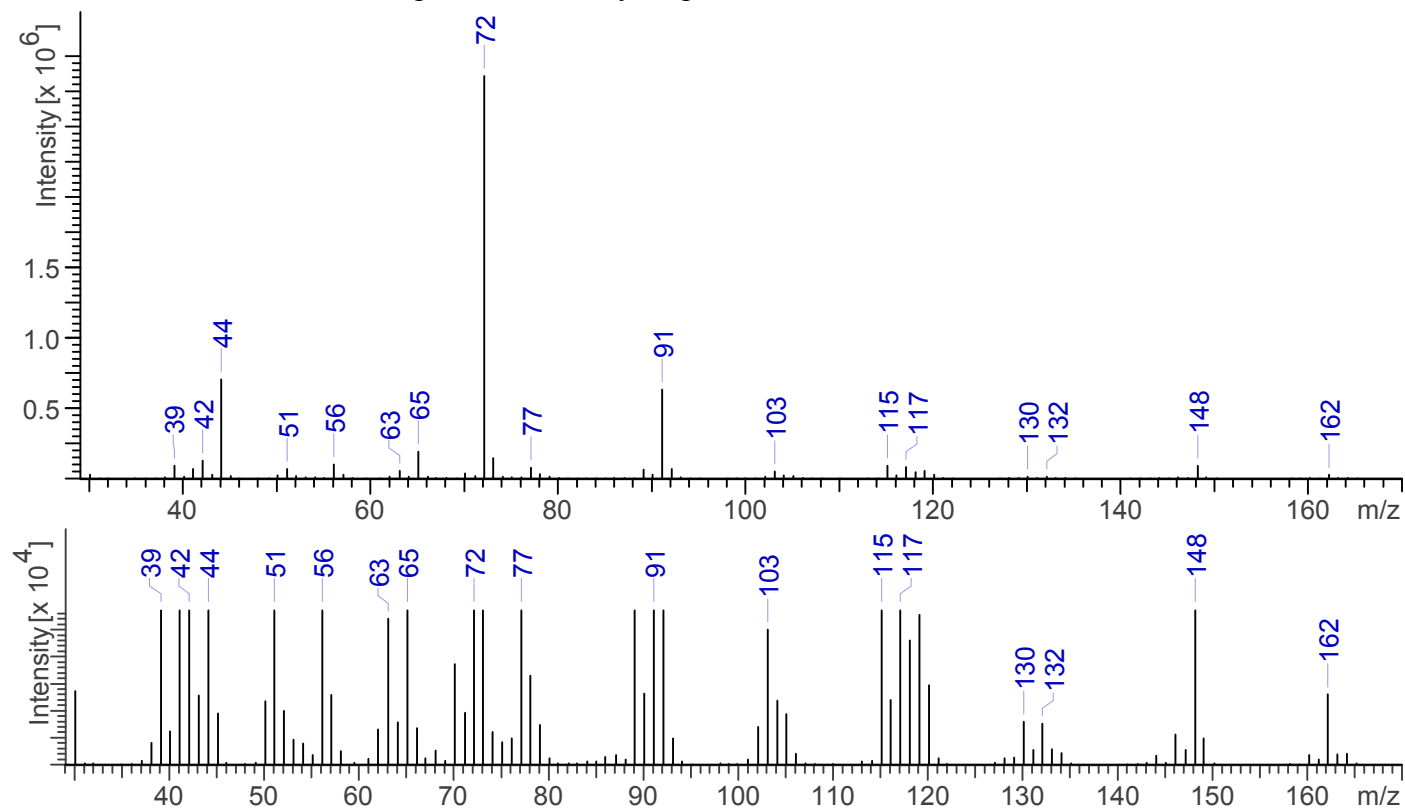

4.3 INFRARED SPECTROSCOPY (FTIR)

Instrument: FTIR with diamond ATR attachment (3 bounce)

Scan Parameters: Number of scans: 32
Number of background scans: 32
Resolution: 4 cm⁻¹
Sample gain: 8
Aperture: 150

FTIR ATR (Diamond, 3 Bounce): N-ethylamphetamine HCl; Lot 2TDM-164-01

