




para-Methoxy furanyl fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.

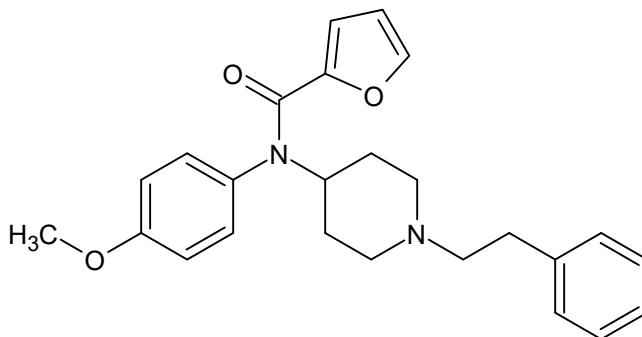

1. GENERAL INFORMATION

IUPAC Name: *N*-(4-methoxyphenyl)-*N*-[1-(phenethylpiperidin-4-yl)furan-2-carboxamide

CAS#: NA

Synonyms: 4-Methoxy Fu-F, *p*-methoxy Fu-F, 4-methoxy furanyl fentanyl, *p*-methoxy furanyl fentanyl, *N*-(4-methoxyphenyl)-*N*-[1-(2-phenylethyl)piperidin-4-yl]furan-2-carboxamide

Source: DEA Reference Material Collection

Appearance: White powder

UV_{max}(nm): Not determined

2. CHEMICAL AND PHYSICAL DATA

2.1 CHEMICAL DATA

Form	Chemical Formula	Molecular Weight	Melting Point (°C)
Base	C ₂₅ H ₂₈ N ₂ O ₃	404.50	Not Determined
HCl	C ₂₅ H ₂₈ N ₂ O ₃ HCl	440.96	Not Determined



para-Methoxy furanyl fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.




3. QUALITATIVE DATA

3.1 NUCLEAR MAGNETIC RESONANCE

Sample Preparation: Dilute analyte to ~14 mg/mL in methanol- d_4 containing TMS for 0 ppm reference and 1,4-BTMSB- d_4 as quantitative internal standard.

Instrument: 400 MHz NMR spectrometer
Parameters: Spectral width: at least containing -3 ppm through 13 ppm
Pulse angle: 90°
Delay between pulses: 45 seconds

^1H NMR: para-Methoxy furanyl fentanyl HCl; Lot# 0546725-1; methanol- d_4 ; 400MHz

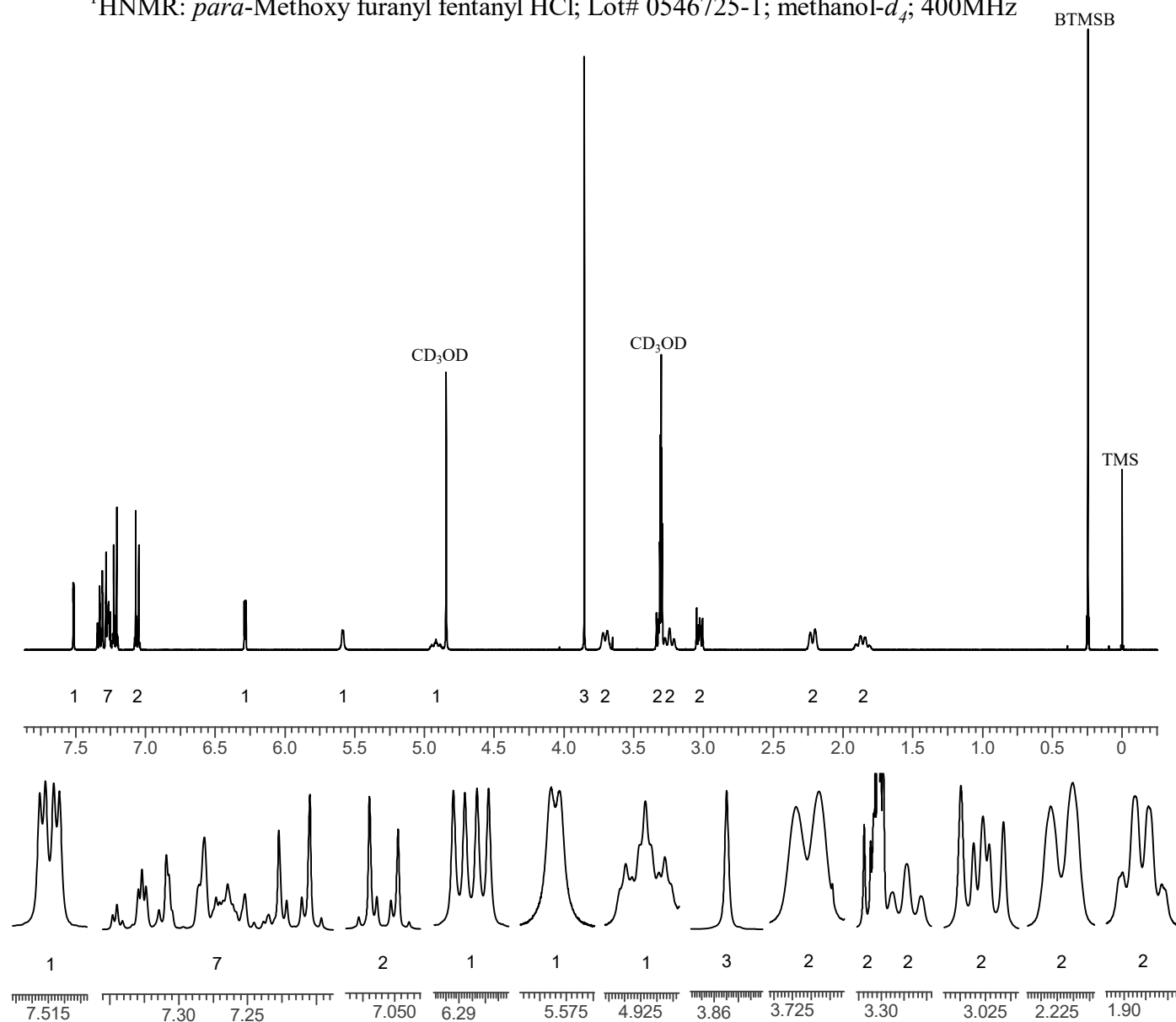



para-Methoxy furanyl fentanyl

The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.



3.2 GAS CHROMATOGRAPHY/MASS SPECTROMETRY

Sample Preparation: Dilute analyte ~4 mg/mL in MeOH

Instrument: Agilent gas chromatograph operated in split mode with MS detector

Column: HP-5 MS (or equivalent); 30m x 0.25 mm x 0.25 μ m

Carrier Gas: Helium at 1.5 mL/min

Temperatures: Injector: 280°C MSD transfer line: 280°C

MS Source: 230°C MS Quad: 150°C

Oven program:

1) 100°C initial temperature for 1.0 min

2) Ramp to 280°C at 12 °C/min

3) Hold final temperature for 9.0 min



Injection Parameters: Split Ratio = 25:1, 1 μ L injected

MS Parameters: Mass scan range: 30-550 amu Threshold: 250

Tune file: stune.u Acquisition mode: scan

Retention Time: 23.44 min

EI Mass Spectrum: *para*-Methoxy furanyl fentanyl HCl; Lot# 0546725-1

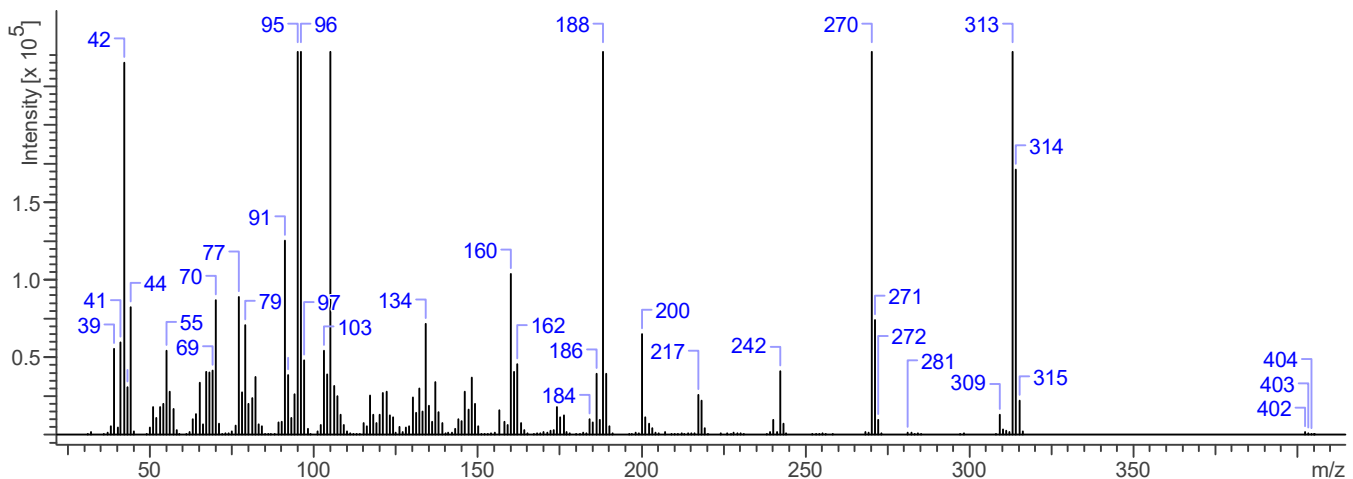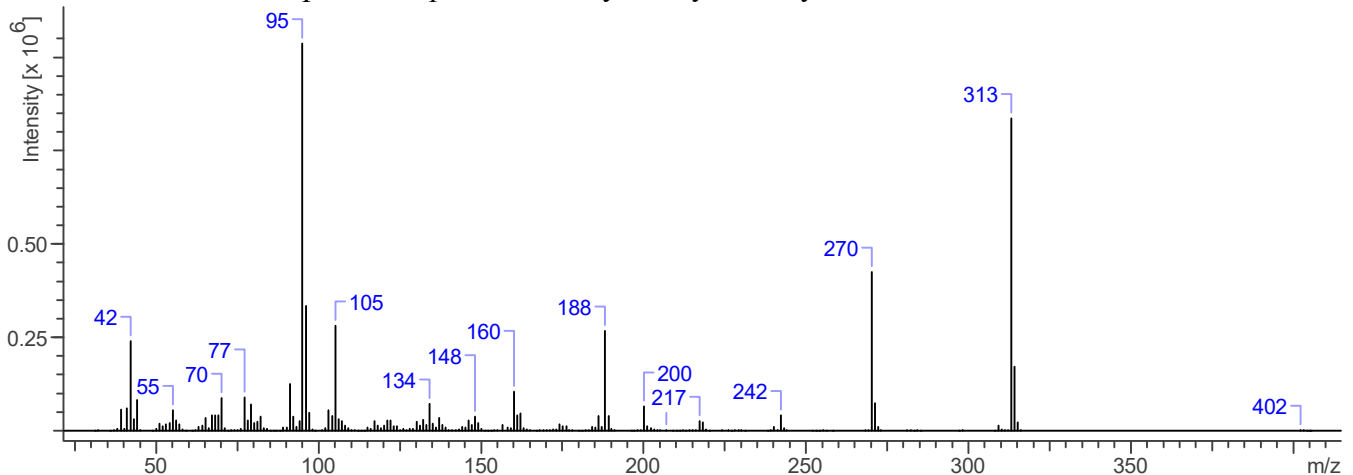



para-Methoxy furanyl fentanyl



The Drug Enforcement Administration's Special Testing and Research Laboratory generated this monograph using structurally confirmed reference material.

3.3 INFRARED SPECTROSCOPY (FTIR)

Instrument: FTIR with diamond ATR attachment (1 bounce)
Scan Parameters: Number of scans: 32
Number of background scans: 32
Resolution: 4 cm⁻¹
Sample gain: 8
Aperture: 80

FTIR ATR (Diamond 1 Bounce): *para*-Methoxy furanyl fentanyl HCl; Lot# 0546725-1

